


4 Gallon 1 Door, 1 Shelf, Self Close, Bench Top Metal Corrosive Safety Cabinet, Blue - CRA1903X

Product Images

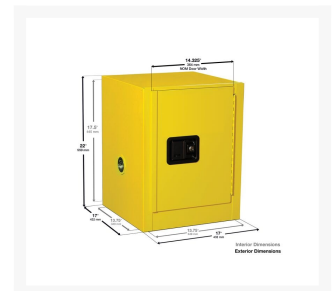

Short Description

- Reinforced double wall construction with 1-1/2 in (3.8 cm) air space on sides, top, bottom, back and doors
- Constructed of 18-gauge steel
- Lead-free, high-gloss, baked-on powder finish
- Continuous piano hinges

- Ground connection
- Four adjustable galvanized leveling legs
- Trilingual warning labels
- 2 in (5.1 cm) raised leakproof sill
- Patented adjustable galvanized steel shelves direct spills to back and bottom of leakproof sump and hold 350 lbs. (158.8 kg) Uniform Distributed Load (UDL)
- Internal and external welded construction
- Two 2 in (5.1 cm) internal vents with fire baffle and cap
- Flush mounted door handle with key lock
- Three-point stainless steel bullet latching system
- 3 in (7.6 cm) centers for shelf brackets
- Self-latching doors
- Door closings: Self-Closing (135° opening)
- Concealed self-close mechanism provides obstruction-free access to top shelf space

Description

Eagle 4 Gallon, self-closing, Bench Top Acid Corrosive Safety Cabinet is constructed of 18-gauge (1 mm) steel. The sides, top, bottom, and doors are double-walled with a 1-1/2 in (3.8 cm) air space between walls. Includes two vents, each have 2 in (5.1 cm) threaded fittings and a fire baffle and cap. Cabinets have a corrosion-resistant powder coat finish and includes trilingual warning label, plus grounding attachment, 3-point latch system, and 2-in (5.1 cm) raised, leakproof door sill. Each shelf is adjustable on 3 in (7.6 cm) centers, supported by 4 brackets, with a 350 lbs (159 kg) safe allowable load. Includes poly shelf and sump trays to help contain small spills and provide extra corrosion protection. Not recommended for use with Sulfuric, Nitric or Hydrochloric acids. Cabinets meet NFPA Code 30 and OSHA requirements.

Specifications

Model No	CRA1903X
Color	Blue
Material Specifications	18-Gauge Steel
FM Approval	Yes
OSHA Compliance	Yes
Complies with NFPA Code 30	Yes
Door Style	Self-close
UPC	48441690121
Dimensions, Exterior	17W x 17Dx 22H
Dimensions, Interior	13.75" W x 13.8" D x 17.5" H (34.9 x 34.9 x 44.5 cm)
Number of Doors	1
Number of Shelves	1
Uniform Fire Code	Yes
Lock	Yes
Gallon Capacity	4
Liter Capacity	15.14
Replacement Shelf Model No.	29935
Pallet Quantity	1
Shelf Depth	13" (33 cm)
Net Weight, lbs	59
Net Weight, kgs	26.76
Adjustable Shelf	Yes

Prop 65

⚠️WARNING: This product can expose you to chemicals including Lead, which is known to the State of

California to cause cancer and birth defects or other reproductive harm. For more information go to www.P65Warnings.ca.gov.