


45 Gallon, 2 Shelves, 2 Door, Sliding Self Close, Metal Corrosive Safety Cabinet, Blue - CRA45X

Product Images

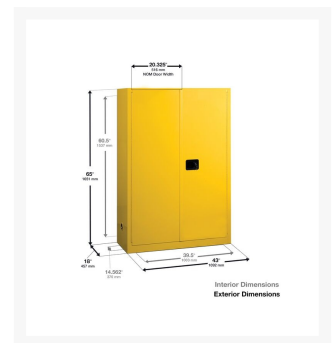

Short Description

- Reinforced double wall construction with 1-1/2 in (3.8 cm) air space on sides, top, bottom, back and doors
- Constructed of 18-gauge steel
- Lead-free, high-gloss, baked-on powder finish
- Continuous piano hinges
- Ground connection

- Four adjustable galvanized leveling legs
- Trilingual warning labels
- 2 in (5.1 cm) raised leakproof sill
- Patented adjustable galvanized steel shelves direct spills to back and bottom of leakproof sump and hold 350 lbs. (158.8 kg) Uniform Distributed Load (UDL)
- Internal and external welded construction
- Two 2 in (5.1 cm) internal vents with fire baffle and cap
- Flush mounted door handle with key lock
- Three-point stainless steel bullet latching system
- 3 in (7.6 cm) centers for shelf brackets
- Self-latching doors
- Door closings: Self-Closing (90° opening)
- Concealed self-close mechanism provides obstruction-free access to top shelf space
- Self-indexing doors guarantee they will close in sequence and ensure a tight closure from top to bottom

Description

Eagle 45 Gallon, sliding self-close, Acid Corrosive Safety Cabinet is constructed of 18-gauge (1 mm) steel. The sides, top, bottom, and doors are double-walled with a 1-1/2 in (3.8 cm) air space between walls. Includes two vents, each have 2 in (5.1 cm) threaded fittings and a fire baffle and cap. Cabinets have a corrosion-resistant powder coat finish and includes trilingual warning label, plus grounding attachment, 3-point latch system, and 2-in (5.1 cm) raised, leakproof door sill. Each shelf is adjustable on 3 in (7.6 cm) centers, supported by 4 brackets, with a 350 lbs (159 kg) safe allowable load. Includes poly shelf and sump trays to help contain small spills and provide extra corrosion protection. Not recommended for use with Sulfuric, Nitric or Hydrochloric acids. Cabinets are FM-approved and meet NFPA Code 30 and OSHA requirements.

Specifications

Model No	CRA45X
Color	Blue
Material Specifications	18-Gauge Steel
FM Approval	Yes
OSHA Compliance	Yes
Complies with NFPA Code 30	Yes
Door Style	Bi-Fold
UPC	048441691708
Dimensions, Exterior	43W x 18Dx 65H
Dimensions, Interior	39.5" W x 14.6" D x 60.5" H (100.3 x 37 x 153.7 cm)
Number of Doors	2
Number of Shelves	2
Uniform Fire Code	Yes
Lock	Yes
Gallon Capacity	45
Liter Capacity	170.34
Replacement Shelf Model No.	29937
Pallet Quantity	1
Shelf Depth	14" (35.6 cm)
Net Weight, lbs	305
Net Weight, kgs	138.35
Adjustable Shelf	Yes

Prop 65

⚠️WARNING: This product can expose you to chemicals including Lead, which is known to the State of

California to cause cancer and birth defects or other reproductive harm. For more information go to www.P65Warnings.ca.gov.